



Low cost, compact, two-channel NMR spectrometer for CP-DNP

Albannay, Mohammed; Johansen, Daniel Højrup; Møllersøe Vinther, Joachim; Zhurbenko, Vitaliy; Ardenkjær-Larsen, Jan Henrik

Publication date:
2017

Document Version
Publisher's PDF, also known as Version of record

[Link back to DTU Orbit](#)

Citation (APA):
Albannay, M., Johansen, D. H., Møllersøe Vinther, J., Zhurbenko, V., & Ardenkjær-Larsen, J. H. (2017). *Low cost, compact, two-channel NMR spectrometer for CP-DNP*. Abstract from EUROMAR 2017, Warsaw, Poland.

General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

Low cost, compact, two-channel NMR spectrometer for CP-DNP

Mohammed M. Albannay¹, Daniel Højrup Johansen¹, Joachim Møllesø Vinther¹, Vitaliy Zhurbenko¹, Jan Henrik Ardenkjaer-Larsen^{1,2}

¹Center for Hyperpolarization in Magnetic Resonance, Department of Electrical Engineering, Technical University of Denmark, 2800 Kgs Lyngby, Denmark. ²GE Healthcare, Brøndby 2605, Denmark.

Hyperpolarization by dissolution dynamic nuclear polarization (dDNP) has become a leading method to boost sensitivity of low gamma nuclei (e.g. ^{13}C , ^6Li , ^{15}N , etc.) in magnetic resonance experiments, typically up to four orders of magnitude¹. Nevertheless, the microwave driven polarization transfer from electron spins to the aforementioned nuclear spins is a slow process, especially at high field (6.7 T)². It has recently been demonstrated that cross-polarization (CP) from protons to low gamma nuclei is beneficial to improve both the build-up time and maximum achievable polarization of the later³. Monitoring of polarization build-up in the solid state at lower frequencies can be achieved with simple low power, single channel NMR spectrometers with low cost. To the best of our knowledge, CP has so far been implemented using dedicated high-resolution spectrometers. We herein present a 10-450 MHz, two-channel, compact, low-cost solution based on a benchtop spectrometer (Kea²; Magritek, New Zealand), a power amplifier (TwinPulse 400, Tomco, Australia) and a home-built duplexer. The improved duplexer allows for experimental flexibility and faster switching between transmitting and receiving modes. All valuable features for hardware demanding NMR sequences, such as CP.

Fabricated on an FR-4 Euroboard PCB, the duplexer installs directly into the Kea² console (Figure 1a) and operates between 10 to 450 MHz. The module accommodates three independently controlled TTL outputs (for power amplifier blanking, microwave gating³ and a spare). The duplexer switching mechanism consists of a $\lambda/4$ cable, PIN diodes (MAP504-1072T; MACOM, USA) and three pairs of crossed 1N4148 diodes thus reducing the dead time $< 2\mu\text{s}$. The preamplifier consists of a single stage low noise amplifier (LNA). Two candidate LNAs were tested (Figure 1b and 1c), AU-2A-150 (Miteq, USA) and a HYPERMAG-developed GaAs LNA and benchmarked against the performance of an INOVA preamplifier (Part No. 01-9104-08; Varian Inc, USA). The selected LNAs provide a maximally flat response with a noise figure < 2 dB. Low gain in the HYPERMAG LNA permits the cascade of the Kea² variable RX amplifier module without saturating the spectrometer's analog-to-digital converter.

The NMR spectrometer was tested at 400 MHz on a liquid sample against an Agilent Direct Drive console. Similar performance was obtained. Likewise, similar performance to standard NMR consoles for solid state NMR detection was obtained. The two channels allow simultaneous pulsing on two nuclei and may therefore be suitable for certain advanced sequences that does not require large memory or fast update rates. A $^1\text{H}/^{13}\text{C}$ CP sequence was implemented and tested with similar performance as an Agilent Direct Drive console.

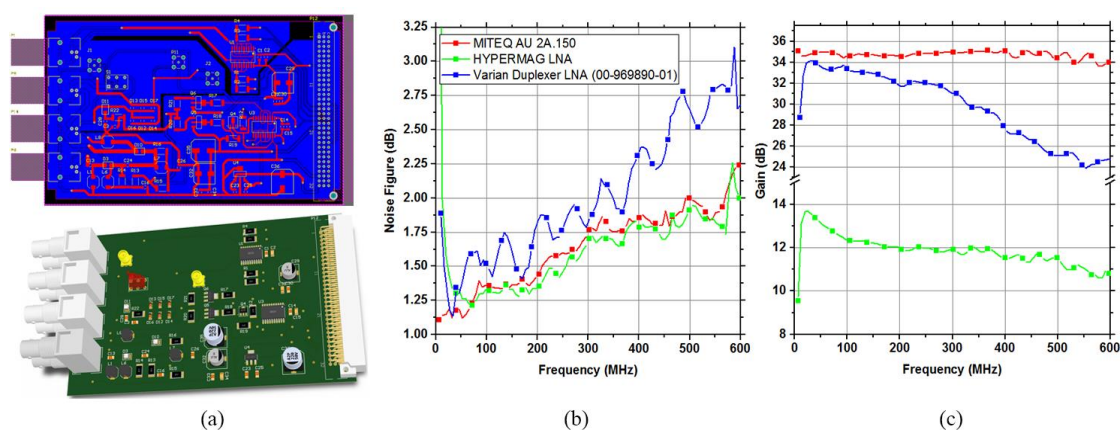


Figure 1(a): Duplexer layout and module (160 mm x 100 mm) with 8 BNC connectors (XMETER (IN), LNA (OUT), Probe (IN/OUT), $\lambda/4$ (IN/OUT), 3 TTL (OUT)). **(1b, 1c):** Measured noise figure and gain of MITEQ AU-2A-150, HYPERMAG and INOVA duplexer LNAs using E4440A spectrometer calibrated with 346B noise source (Keysight Technologies, USA).

References

1. Ardenkjaer-Larsen, J. H. *et al.* Increase in signal-to-noise ratio of $> 10,000$ times in liquid-state NMR. *Proc. Natl. Acad. Sci. U. S. A.* **100**, 10158–10163 (2003).
2. Yoshihara, H. A. I. *et al.* High-field dissolution dynamic nuclear polarization of [1- ^{13}C]pyruvic acid. *Phys. Chem. Chem. Phys.* **18**, 12409–12413 (2016).
3. Bornet, A. *et al.* Microwave-gated dynamic nuclear polarization. *Phys. Chem. Chem. Phys.* **18**, 30530–30535 (2016).